

# AIV Ab ELISA

## Validated from

- \*IZS delle Venezie, Legnaro (Padova), Italy
- \*FGI «ARRIAH», Federal Center For Animal Health, Russia
- \*University of Hokkaido, Graduated School of Veterinary Medicine

## Avian Influenza type A virus antibody

Wild birds and some poultry infected by HPAI show no clinical signs, and their virus shedding amount is too small to be detected by an antigen capture immunoassay. The AniGen AIV Ab ELISA has been developed for screening these latent AIV carriers by detecting AIV antibodies.



## Indications


- Avian influenza virus antibody detection
- Monitoring of Avian Influenza virus carrier
- Antibody detection for migratory bird surveillance

## Special Features

- No cross-reaction with other avian virus positive sera
- Applicable to various species
- Specimen: Serum, plasma or egg yolk
- No sample dilution is required
- Reading time: 45 minutes
- Sensitivity: Chicken 98.2%, Duck 97.5%, Turkey 97.1%
- Specificity: Chicken 97.3%, Duck 93.8%, Turkey 100%

## Quick Procedure

1. Prepare AIV NP antigen coated test plate
2. Add 50  $\mu$ l of controls and sample to wells
3. Add 50  $\mu$ l of anti AIV antibody-HRP to each wells
4. Incubate plate for 30 minutes at 37  $^{\circ}$ C
5. Wash plate 6 times
6. Add 100  $\mu$ l of substrate (Ready to use) and incubated the wells for 10 minutes at room temperature
7. Add 100  $\mu$ l of stopping solution
8. Measure the optical density (OD) at 450 nm with reference wavelength at 620nm
9. PI value=[1-(OD sample/mean OD negative)] x 100

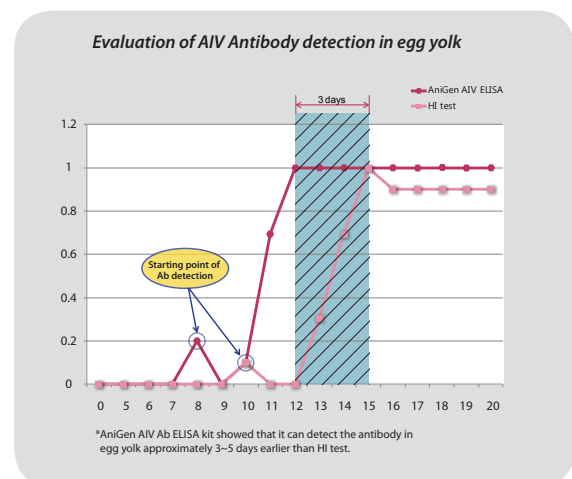

## Ordering Information

Cat. No.	Description	Type	Pack size
EB45-02	AIV Ab ELISA	Microplate	480 Wells/Kit