

FCoV Ab

General Information

Purpose : Detection of Feline Corona virus antibody

Principle : Immunochromatographic assay

Specimen : Whole blood, plasma or serum

Materials provided :

- Anigen Rapid FCoV Ab test device
- Assay diluents bottle
- Anticoagulant tube
- Disposable capillary tube

Special Features :


- One-step testing procedure : Quick & Accurate
- High sensitivity and specificity
 - Sensitivity : 96.0 % vs. IFA
 - Specificity : 97.9 % vs. IFA



◆ **How to do FIP diagnosis with Rapid FCoV Ab test**

- Many healthy cats and cats with conditions other than FIP are FCoV Ab positive.
- Some cats with effusive FIP appear to have low titers or are FCoV Ab negative.
- FCoV Ab testing should only be performed to determine FIP in conjunction with a compatible history, clinical signs, and examination of effusions or blood for high globulins and a low A:G ratio.

Test Procedure

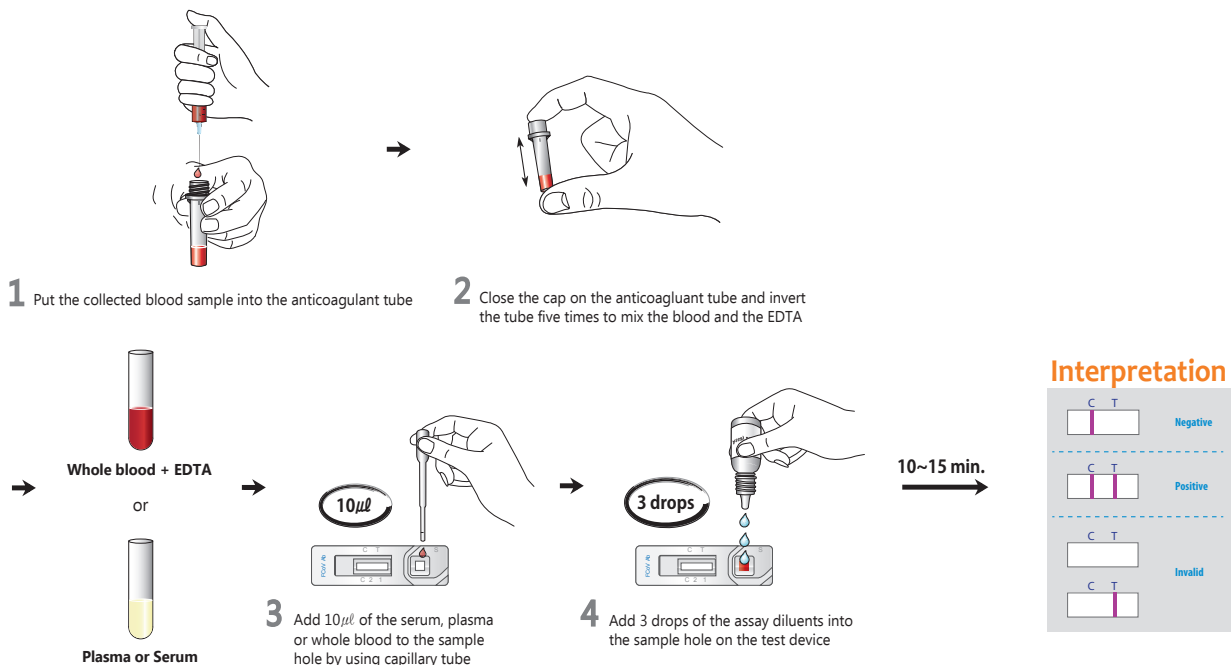

Ordering Information

Cat. No.	Description	Type	Pack size
RB2203DD	FCoV Ab	Device	10 Tests/Kit